
# beurer

## cellulite releaZer®
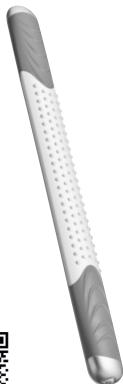
|           |                                                     |     |
|-----------|-----------------------------------------------------|-----|
| <b>DE</b> | cellulite releaZer<br>Gebrauchsanweisung.....       | 2   |
| <b>EN</b> | cellulite releaZer<br>Instructions for use.....     | 20  |
| <b>FR</b> | cellulite releaZer<br>Mode d'emploi .....           | 39  |
| <b>ES</b> | cellulite releaZer<br>Instrucciones de uso .....    | 58  |
| <b>IT</b> | cellulite releaZer<br>Istruzioni per l'uso .....    | 77  |
| <b>TR</b> | cellulite releaZer<br>Kullanım kılavuzu .....       | 96  |
| <b>RU</b> | cellulite releaZer<br>Инструкция по применению .... | 115 |
| <b>PL</b> | cellulite releaZer<br>Instrukcja obsługi .....      | 134 |

# DEUTSCH

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# ENGLISH

**Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.**

## **Dear customer,**

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, gentle therapy, blood pressure/diagnosis, weight, massage, beauty, air and baby care.

With kind regards,  
Your Beurer team.

## **Table of contents**

|                                                                                  |           |
|----------------------------------------------------------------------------------|-----------|
| <b>1. Getting to know your device .....</b>                                      | <b>23</b> |
| <b>2. Signs and symbols .....</b>                                                | <b>24</b> |
| <b>3. Intended use .....</b>                                                     | <b>25</b> |
| <b>4. Notes .....</b>                                                            | <b>26</b> |
| <b>5. Device description .....</b>                                               | <b>28</b> |
| <b>6. Initial use/charging process .....</b>                                     | <b>29</b> |
| <b>7. Use .....</b>                                                              | <b>30</b> |
| 7.1 Preparing for use .....                                                      | 30        |
| 7.2 Switching on the cellulite releaZer / selecting vibrati-<br>on setting ..... | 31        |
| 7.3 Tips for use .....                                                           | 32        |
| 7.4 Treatment on the outside of the thigh .....                                  | 32        |
| 7.5 Treatment on the rear of the thigh .....                                     | 34        |
| <b>8. Cleaning and maintenance .....</b>                                         | <b>35</b> |
| <b>9. Disposal .....</b>                                                         | <b>36</b> |
| <b>10. Technical data .....</b>                                                  | <b>37</b> |
| <b>11. Warranty / service .....</b>                                              | <b>37</b> |

## **WARNING**

- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

## **Included in delivery**

- 1x cellulite releaZer / 1x micro-USB charging cable / 1x mains adapter / 1x these instructions for use

# **1. Getting to know your device**

## **What is cellulite?**

Cellulite, also known as orange peel skin, is something that lots of women have. Oestrogen, the female sex hormone is to blame. It is primarily responsible for orange peel skin with its dimpling and bumps on the bottom and thighs. Cellulite is not an illness and it normally isn't painful. Cellulite is a purely cosmetic problem.

## **What causes cellulite?**

Humans have lobules in their subcutaneous skin and these lobules are separated from one another by septae.

In women, the septae run perpendicular to the surface of the skin. This means that women have softer and more flexible skin than men. The disadvantage of this is that women's softer skin provides less resistance to the pressure of the fatty tissue below. The subcutaneous fat therefore penetrates into the top layers of skin more easily and makes the skin bulge out visibly.

## **How does the cellulite releaZer work?**

When used correctly, the cellulite releaZer increases microcirculation in the tissue by using vibration.

## 2. Signs and symbols

The following symbols appear in these instructions for use:

|                                                                                  |                                                                                               |
|----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------|
|  | <b>WARNING</b><br>Warning notice indicating a risk of injury or damage to health              |
|  | <b>IMPORTANT</b><br>Safety note indicating possible damage to the device/accessory            |
|  | <b>Note</b><br>Note on important information                                                  |
|  | Observe the instructions for use                                                              |
|  | Manufacturer                                                                                  |
|  | Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE |
|  | Only use indoors                                                                              |

The device is double protected and therefore corresponds to protection class 2

### 3. Intended use

#### **WARNING**

Only use the cellulite releaZer for massage treatment of cellulite. The cellulite releaZer is intended for private use only. The cellulite releaZer is not a medical device; rather it is a massage treatment device. The cellulite releaZer may only be used for the purpose for which it was designed and in the manner specified in the instructions for use. Any form of improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or incorrect use.

- Do not use the cellulite releaZer if you suffer from a medical abnormality or injury in the area to be massaged (e.g. open wounds).
- Do not use the cellulite releaZer during pregnancy.
- Do not use the cellulite releaZer whilst sleeping.
- Do not use the cellulite releaZer in a vehicle.
- Do not use the cellulite releaZer on animals.
- Do not use the cellulite releaZer whilst undertaking any activity where an unexpected reaction could be dangerous.
- Do not use the cellulite releaZer on your head.
- Do not use the cellulite releaZer on your face.
- Do not use the cellulite releaZer on your genitals.
- Do not use the cellulite releaZer after the consumption of substances that lead to impaired cognitive function (e.g. pain medication, alcohol).

- Do not use the cellulite releaZer for longer than 15 minutes (risk of overheating) and leave the device to cool down for at least 15 minutes before using it again. The cellulite releaZer has an automatic switch-off function meaning that it switches off automatically after 15 minutes of use.

If you are unsure whether the cellulite releaZer is suitable for you, please consult your doctor.

Please consult your doctor before using the cellulite releaZer

- if you suffer from a severe illness or have recently had surgery,
- if you have a pacemaker, implant or other aids,
- in the case of feverish colds, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nerve disorders (e.g. sciatica) or acute inflammations,
- if you are from suffering pain, the cause of which is unknown.

## 4. Notes

### WARNING

#### **Safety notes regarding the use of the cellulite releaZer**

- Keep the packaging material away from children. Choking hazard!
- The cellulite releaZer is not a medical device. It is a massage treatment device.
- Before use, ensure that the device and accessories have no visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.
- Repairs must only be carried out by Customer Services or authorised retailers. Under no circumstances should you open or repair the cellulite releaZer yourself, as faultless functionality can no longer be guaranteed thereafter.


Failure to comply with this instruction will void the warranty. Never attempt to repair the device yourself!

- Only use the device with the supplied accessories.

## **Notes on handling batteries**

- If your skin or eyes come into contact with fluid from a rechargeable battery cell, flush out the affected areas with water and seek medical assistance.
- ⚠ Choking hazard! Small children may swallow and choke on batteries. Store batteries out of the reach of small children.
- Observe the plus (+) and minus (-) polarity signs.
- Protect batteries from excessive heat.
- ⚠ Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, split or crush the batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter 6. Initial use/charging process).
- In order to achieve as long a battery service life as possible, fully charge the battery at least twice per year.

## 5. Device description

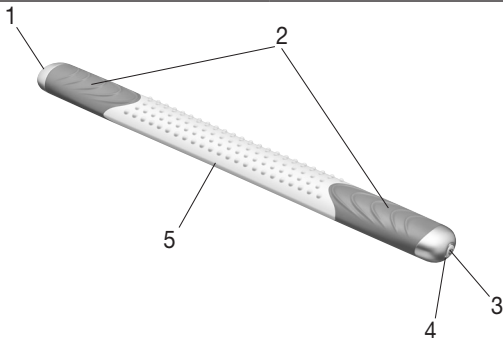1. Micro USB connection

2. Handles

3. ON/OFF button

4. Function light

5. Four treatment surfaces  
(smooth, nubs, blunt edge, and  
sharp edge)

| Function light status  | Charging function<br>(charging cable plugged in) |
|--------------------------------------------------------------------------------------------------------|--------------------------------------------------|
| Function lights flash in a circle                                                                      | Charging                                         |
| Function lights illuminate continuously                                                                | Charging complete                                |
| The left and right function lights flash alternately                                                   | Low battery                                      |

## 6. Initial use/charging process

### Note

Before using the cellulite releaZer for the first time, let it charge for a minimum of 3.5 hours. Proceed as follows:

- Connect the charging cable to the mains adapter and the cellulite releaZer (*Figure 1*). Then insert the mains adapter into a suitable socket. Alternatively, you can also charge the cellulite releaZer with the charging cable via the USB port on your computer or laptop (*Figure 2*).

When the device is charging, the function lights flash in a circle. When the cellulite releaZer is fully charged, the function lights illuminate continuously. If the left and right function lights flash alternately during operation, you have to charge the cellulite releaZer.

Figure 1


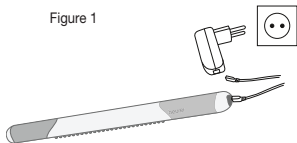
Figure 2
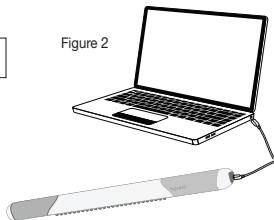
## 7. Use

### 7.1 Preparing for use

Before you start using the device, you should gently rub massage oil onto the parts of the body to be treated.

Do not use any oils or creams with a strong cooling or heating effect, as they may cause skin irritation.

1. Put a couple of drops of massage oil onto the palm of your hand.
2. Now rub the massage oil evenly onto the parts of the body to be massaged.

- Before treatment, wipe your hands with a towel to prevent them from slipping when using the device and to prevent damage to the silicone surface on the handles.

## **7.2 Switching on the cellulite releaZer / selecting vibration setting**

Press the ON/OFF button to switch on the cellulite releaZer. The function light now illuminates blue and the first vibration setting is activated. The cellulite releaZer has a total of four different vibration settings.

| <b>Action</b>           |   | <b>Vibration setting</b>                                                                           |
|-------------------------|---|----------------------------------------------------------------------------------------------------|
| 1st press of the button | = | The cellulite releaZer switches on and the function light illuminates blue. Low vibration setting. |
| 2nd press of the button | = | The second function light illuminates blue. Medium vibration setting.                              |
| 3rd press of the button | = | The third function light illuminates blue. High vibration setting.                                 |
| 4th press of the button | = | The fourth function light flashes blue. Intermittent vibration setting.                            |
| 5th press of the button | = | The cellulite releaZer switches off and the function light goes out.                               |

## **7.3 Tips for use**

- During use, hold the cellulite releaZer with both hands using the handles.
- For optimum results, use the cellulite releaZer every two days.
- Severe reddening or bruises may occur at the treated areas after use. This is completely normal and indicates that you have used the cellulite releaZer correctly.

### **Note**



**You can use the cellulite releaZer on all areas where you have cellulite. The examples in the following chapter showing the outside of the thigh (iliotibial band) and the rear of the thigh (biceps femoris muscle) demonstrate how the cellulite releaZer should be used.**

## **7.4 Treatment on the outside of the thigh**

The treatment time for the outside of the thigh is approx. 7 minutes and consists of a total of three treatment phases. To start treatment, proceed as follows:

1. Follow the instructions in the "Preparing for use" chapter.
2. Place the leg to be treated on a chair, stool or the edge of the bath tub, for example. Carry out the entire treatment on the thigh whilst standing.
3. Select your desired vibration setting (1, 2, 3 or 4) using the ON/OFF button. The cellulite releaZer starts to vibrate.


4. Treatment phase I: Treat the outside of the thigh with the nub side using small, fast, circular movements and moderate pressure (for approx. 3 minutes).
5. Treatment phase II: Treat the outside of the thigh with the nub side using slow, deep, circular movements and increased pressure (for approx. 3 minutes).
6. Treatment phase III: The third treatment phase is drainage. Glide the device towards the centre of the body/heart using the coloured edge and moderate pressure. When doing so, start the movement at the upper part of the outside of the thigh and slowly move the starting point down towards the knee. One stroke should be approximately 10 cm long. Carry out this process 5 times moving down the outside of the thigh until you reach the knee area. At the end, move the device over the tissue of the entire outside of the thigh once.

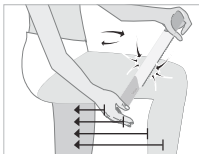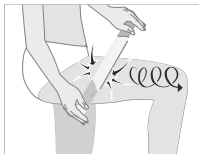
## **7.5 Treatment on the rear of the thigh**

The treatment time for the rear of the thigh and buttocks is approx. 7 minutes and consists of a total of three treatment phases. To start treatment, proceed as follows:

1. Follow the instructions in the "Preparing for use" chapter.
2. Place the leg to be treated on a chair or stool. Carry out the entire treatment on the thigh whilst standing.
3. Select your desired vibration setting (1, 2, 3 or 4) using the ON/OFF button. The cellulite releaZer starts to vibrate.
4. Treatment phase I: Treat the rear of the thigh with the nub side using small, fast, circular movements and moderate pressure (for approx. 3 minutes).
5. Treatment phase II: Treat the rear of the thigh with the nub side using slow, deep, circular movements and increased pressure (for approx. 3 minutes).
6. Treatment phase III: The third treatment phase is drainage. Glide the cellulite releaZer towards the centre of the body/heart using the coloured edge and moderate pressure. When doing so, start the movement at the upper part of the rear of the thigh (buttocks) and slowly move the starting point down towards the knee. One stroke should be approximately 10 cm

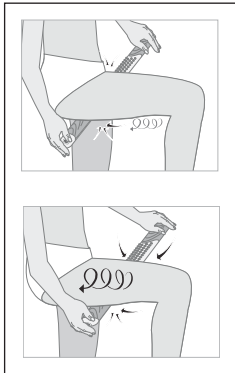long. Carry out this process 5 times moving down the rear of the thigh until you reach the knee area. At the end, move the device over the tissue of the entire rear of the thigh once.
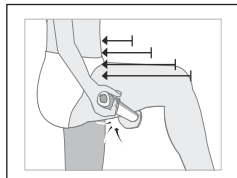
## 8. Cleaning and maintenance

The service life of the cellulite releaZer depends on careful use.

### **IMPORTANT**

- Clean the cellulite releaZer with a wet or dry cloth after each use. Do not use any abrasive cleaning products.
- Protect the cellulite releaZer from knocks, dust, chemicals, drastic changes in temperature, electromagnetic fields and nearby sources of heat (ovens, heaters).

## 9. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life.

Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

### Battery disposal


The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries. To remove the battery from the device, pry open the cover on the micro USB connection using a screwdriver. Loosen the four screws on the battery holder. Then remove the battery.

The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury
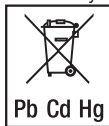
## 10. Technical data

|                                                                    |                                                           |
|--------------------------------------------------------------------|-----------------------------------------------------------|
| Supply                                                             | Input: 100-240V~; 50/60 Hz; 0.2 A<br>Output: 5.0V ---1.0A |
| Dimensions                                                         | 48 x 5 x 3 cm                                             |
| Weight                                                             | Approx. 570 g                                             |
| <u>Battery:</u><br>Capacity<br>Nominal voltage<br>Type designation | 2200 mAh<br>3.7 V<br>Lithium-ion                          |

## 11. Warranty / service

In case of a claim under the warranty please contact your local dealer or the local representation which is mentioned in the list “service international”.  
In case of returning the unit please add a copy of your receipt and a short report of the defect.

The following warranty terms shall apply:

1. The warranty period for BEURER products is either 3 years or –if longer– the country-specific warranty period from date of purchase.  
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.

2. Repairs (complete unit or parts of the unit) do not extend the warranty period.
3. The warranty shall not be valid for damages because of
  - a. improper treatment, e.g. nonobservance of the user instructions.
  - b. repairs or tampering by the customer or unauthorised third parties.
  - c. transport from the manufacturer to the consumer or during transport to the service centre.
  - d. The warranty shall not be valid for accessories which are subject to normal wear and tear (cuff, batteries etc.).
4. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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